


# Morrison

STEERABLE NEEDLE®

Your first choice for  
deep-seated targets

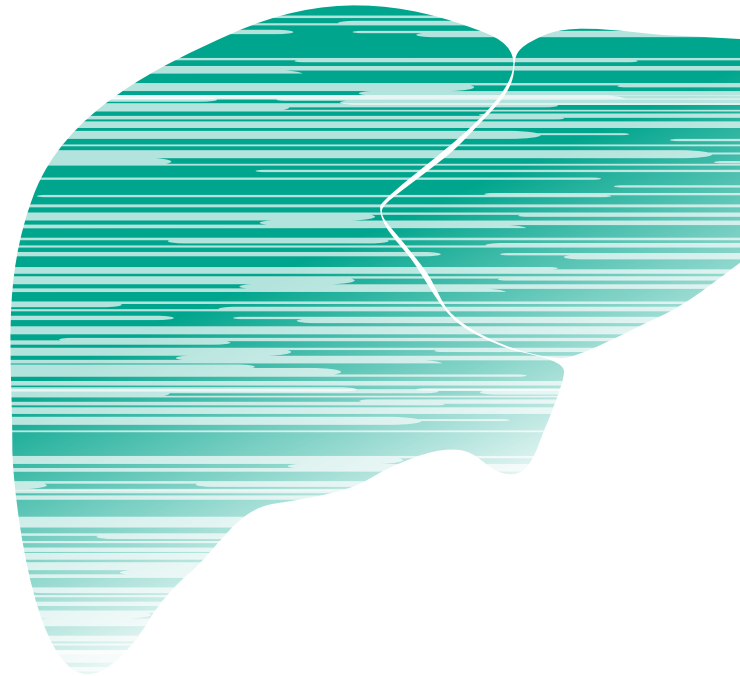

## Indication for use:

Musculoskeletal percutaneous injection and aspiration procedures, and tissue sampling.




**AprioMed**

- Steer distal end of needle with live feedback under image guidance
- Facilitates precise needle placement
- Curve around sensitive structures or other obstacles

**USE MORRISON STEERABLE NEEDLE® TO:**

- Make major adjustments around bones, organs or neurovascular structures
- Make minor adjustments near target
- Make multiple adjustments as needle is advanced

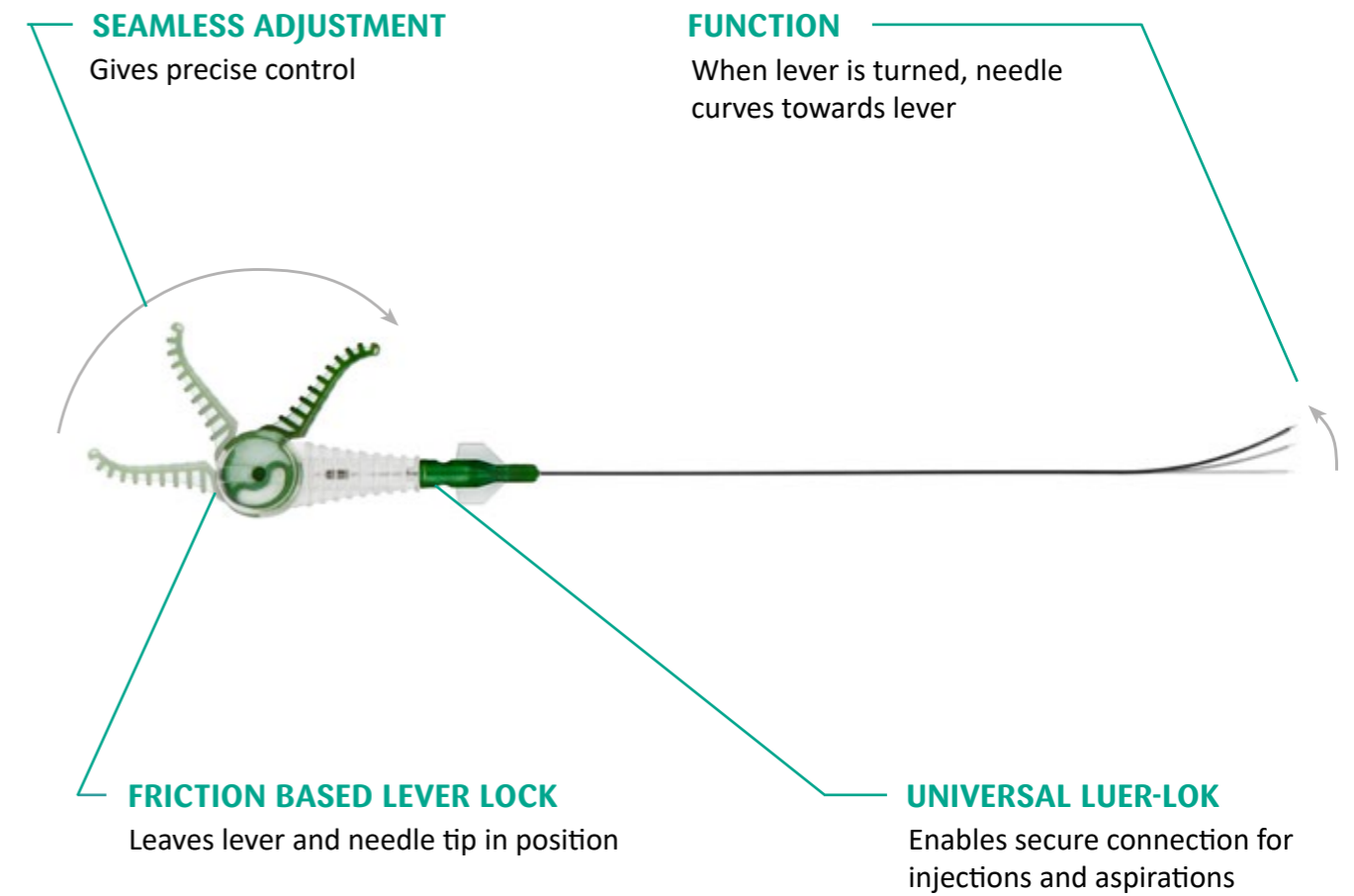


**PENCIL POINT STYLET TIP**

Facilitates straight and smooth needle advancement

**30° BACK BEVEL NEEDLE TIP**

Provides perfect cut through tissue

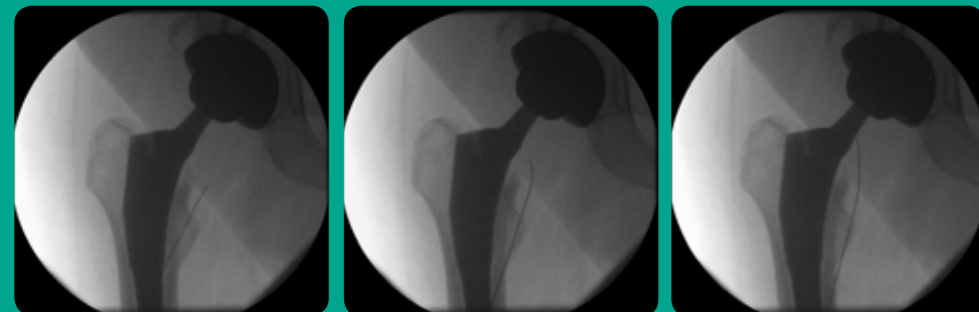

Deep-seated target access


Morrison Steerable Needle™ used for access to a deep-seated target in the hip joint.

## TECHNICAL SPECIFICATIONS

### Morrison Steerable Needle<sup>®</sup>, 21G


Package contents:

- Steerable Stylet
- Needle
  - Inner diameter 0.68 mm
  - Outer diameter 0.82 mm (21G)
  - Length 167 mm

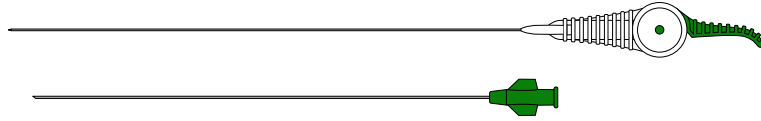


## PERFORMANCE SPECIFICATIONS

- a) Length, upright position
- b) Deflectable tip length
- c) Maximum tip deflection
- d) Curve radius
- e) Seamless adjustment

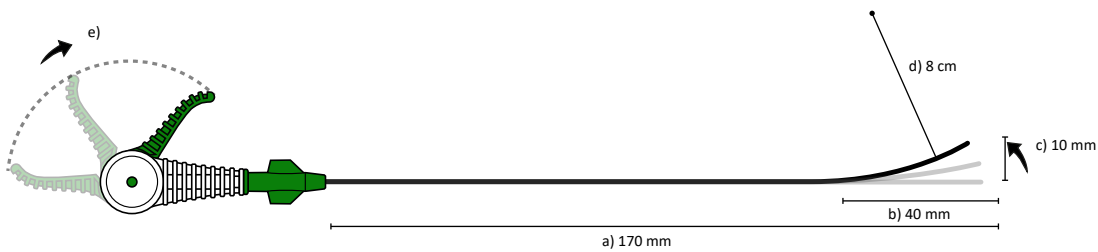

## ORDERING INFORMATION

Product number	Product description	Minimum Order Quantity
13500	Morrison Steerable Needle <sup>®</sup> 21G 17 cm	10 units

 Indication for use: Musculoskeletal percutaneous injection and aspiration procedures, and tissue sampling. Before using Morrison Steerable Needle<sup>®</sup> read the instructions for use which accompany the product for indications, contraindications, warnings and precautions. Morrison Steerable Needle<sup>®</sup> is a trademark of AprioMed AB. Patents pending.

Not available in the USA.

### HEAD OFFICE

AprioMed AB  
Viridings Allé 28  
SE-754 50 Uppsala, Sweden  
Phone: +46 18 430 14 40  
Free Fax (Int'l): +46 770 339972

### US SALES OFFICE

AprioMed, Inc.  
2 Palmer Dr, Suite 1  
Londonderry, NH 03053, USA  
Phone: +1 603 421 0875  
Fax (US toll free): 866 239 1797



info@apriomed.com / www.apriomed.com

APRIOMED develops, manufactures, markets and sells innovative medical devices and related services within the field of interventional radiology. We aim to deliver, in close collaboration with healthcare practitioners, innovative quality tools to achieve the most favourable solutions for radiologists worldwide. AprioMed is headquartered in Uppsala, Sweden, with a sales subsidiary in Londonderry, NH, USA. Our products are sold in more than 40 countries by AprioMed (direct sales in US, UK, Ireland and Sweden) and our local distributors.