
HIGHLIGHTED FEATURES

Gangi-SoftGuard[®]

Coaxial Needle

Tactility, safety and visual feedback...with a spring




AprioMed

DEVICE DESCRIPTION

- Gangi-SoftGuard® combines the well-known technique of a coaxial biopsy needle with an all new spring loaded blunt tip stylet providing extra control and safety during soft tissue biopsy.
-

SPRING BACK PIN

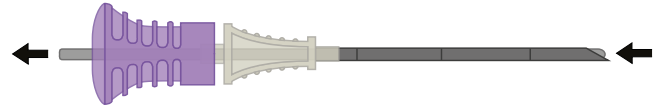

- Pin springs back simultaneously as blunt tip springs back
 - Provides visual and tactile feedback of the blunt tip
 - Pressing the pin gives control of the blunt tip
-



BEVEL INDICATOR



- Marking on hub indicates bevel position
-

BIOPSY NEEDLE

- Provides the ability to cut through tissue
- Centimeter markings on needle for easy depth evaluation

SPRING LOADED BLUNT STYLET

- Enables penetration without removal of the blunt needle
- Provides visual and tactile feedback

BLUNT TIP

- Reduces the risk of inadvertent perforation of adjacent structures when advancing the needle




TECHNICAL SPECIFICATIONS

Package contents:

- Blunt stylet
- Trocar stylet
- Needle
- Reference guide

Gangi-SoftGuard[®], 12G

Needle inner diameter: 2.48 mm
Available needle length: 12.1 cm
Suitable biopsy device size: 13G

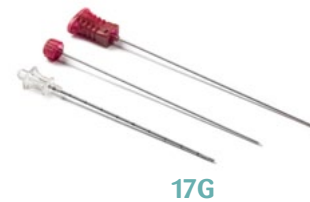

Gangi-SoftGuard[®], 13G

Needle inner diameter: 2.18 mm
Available needle lengths: 7.1 cm, 12.1 cm, 17.1 cm
Suitable biopsy device size: 14G

Gangi-SoftGuard[®], 15G

Needle inner diameter: 1.73 mm
Available needle lengths: 7.1 cm, 12.1 cm, 17.1 cm
Suitable biopsy device size: 16G

Gangi-SoftGuard[®], 17G

Needle inner diameter: 1.33 mm
Available needle lengths: 7.1 cm, 12.1 cm, 17.1 cm
Suitable biopsy device size: 18G

ORDERING INFORMATION

Product number	Product description	Minimum Order Quantity
SG12-121	Gangi-SoftGuard [®] 12G 12.1 cm	5 units
SG13-071	Gangi-SoftGuard [®] 13G 7.1 cm	5 units
SG13-121	Gangi-SoftGuard [®] 13G 12.1 cm	5 units
SG13-171	Gangi-SoftGuard [®] 13G 17.1 cm	5 units
SG15-071	Gangi-SoftGuard [®] 15G 7.1 cm	5 units
SG15-121	Gangi-SoftGuard [®] 15G 12.1 cm	5 units
SG15-171	Gangi-SoftGuard [®] 15G 17.1 cm	5 units
SG17-071	Gangi-SoftGuard [®] 17G 7.1 cm	5 units
SG17-121	Gangi-SoftGuard [®] 17G 12.1 cm	5 units
SG17-171	Gangi-SoftGuard [®] 17G 17.1 cm	5 units



Indication for use: As a guiding needle in obtaining core biopsy samples from soft tissue such as liver, kidney, spleen, lymph nodes and various soft tissue lesions. Before using Gangi-SoftGuard[®] read the instructions for use which accompany the product for indications, contraindications, warnings and precautions. Gangi-SoftGuard[®] is a registered trademark of AprioMed AB. Patents pending.

29375 Rev 01 2016-06

HEAD OFFICE
AprioMed AB
Virdings Allé 28
SE-754 50 Uppsala, Sweden
Phone: +46 18 430 14 40
Free Fax (Int'l): +46 770 339972

E-mail: info@apriomed.com
www.apriomed.com

AprioMed

US SALES OFFICE
AprioMed, Inc.
2 Palmer Dr, Suite 1
Londonderry, NH 03053, USA
Phone: +1 603 421 0875
Fax (US toll free): 866 239 1797